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
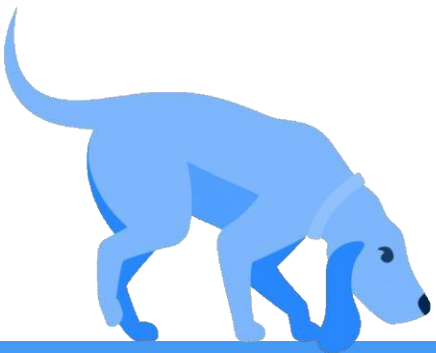
aleatorizado

Agregar evento

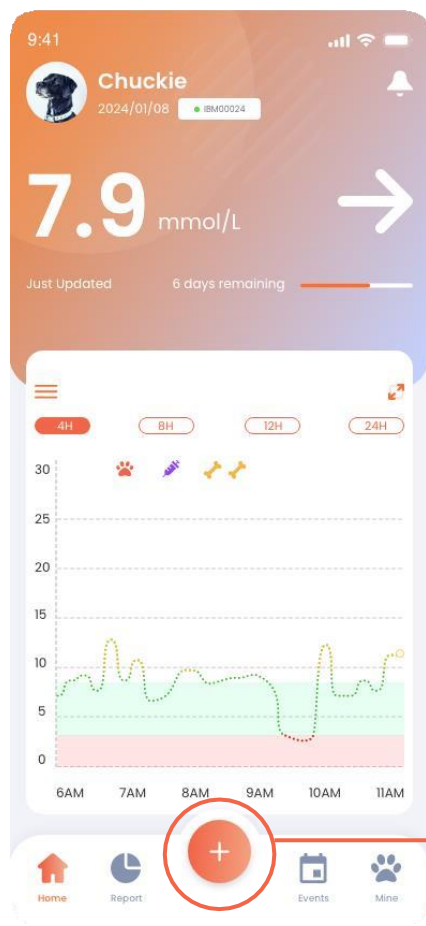
This form is titled "Add new event" and has a close button (X). It features five icons: Exercise, Insulin, Food, BG, and Weight. The "BG" icon is highlighted. Below the icons, there is a section for "Please Enter BG Value" with a text input field. A "Photo Record" section includes a camera icon. At the bottom, there is a date and time selector showing "2021-06-01 14:08" and a "Confirm" button.This form is titled "Add new event" and has a close button (X). It features five icons: Exercise, Insulin, Food, BG, and Weight. The "Food" icon is highlighted. Below the icons, there is a section for "Food Name" with a text input field. An "Amount" section includes a text input field and a unit selector (g). At the bottom, there is a date and time selector showing "2021-06-01 14:08" and a "Confirm" button.This form is titled "Add new event" and has a close button (X). It features five icons: Exercise, Insulin, Food, BG, and Weight. The "Exercise" icon is highlighted. Below the icons, there is a section for "Please enter" with a text input field. A "Remark" section includes a text input field. At the bottom, there is a date and time selector showing "2021-06-01 14:08" and a "Confirm" button.This form is titled "Add new event" and has a close button (X). It features five icons: Exercise, Insulin, Food, BG, and Weight. The "Insulin" icon is highlighted. Below the icons, there is a section for "Insulin Type" with a dropdown menu. An "Insulin Dose" section includes a text input field and a unit selector (Units). At the bottom, there is a date and time selector showing "2021-06-01 14:08" and a "Confirm" button.

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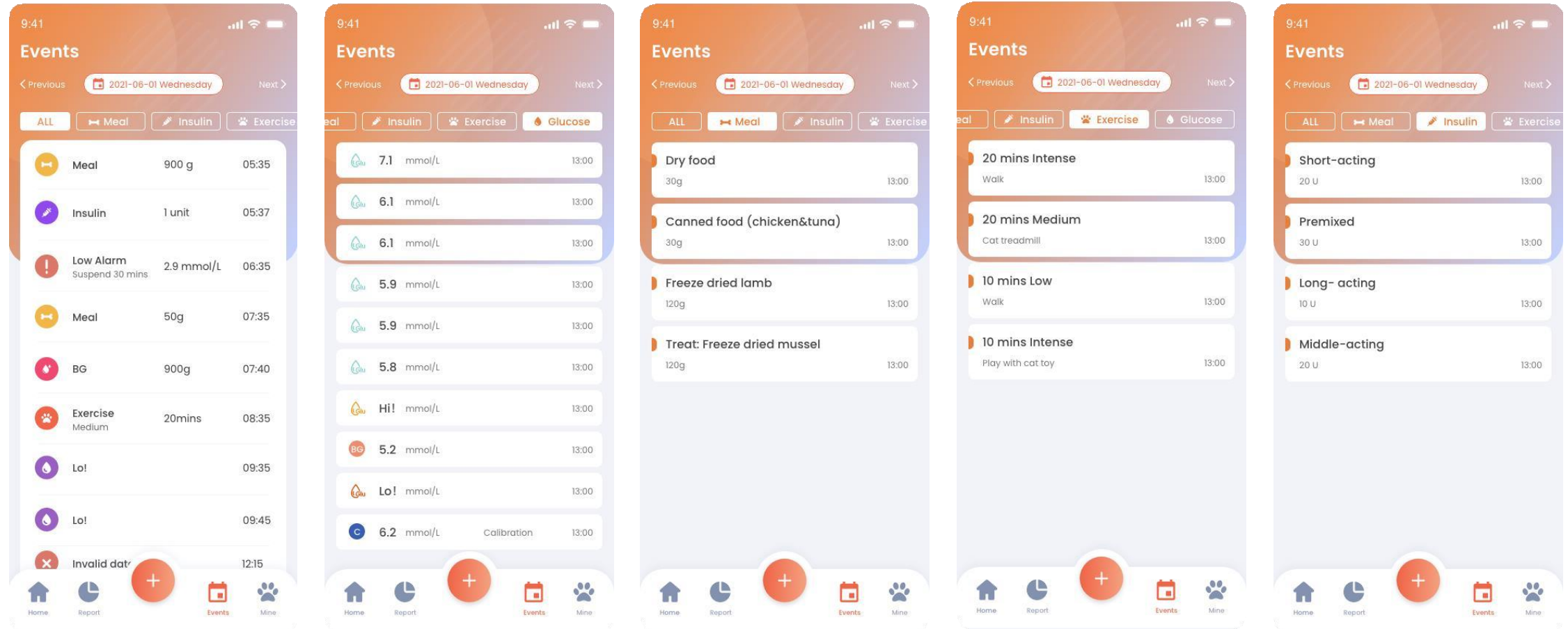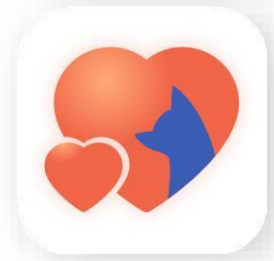
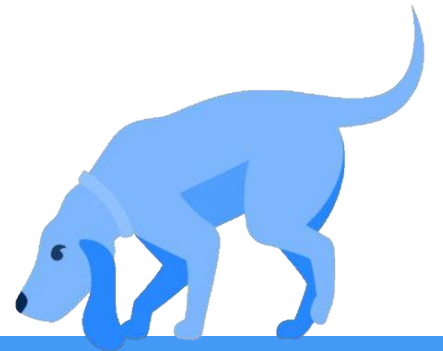
A screenshot of a mobile phone screen displaying the FurLink app's sign-up interface. The status bar at the top shows the time 17:22, signal strength, Wi-Fi, and battery level at 61%. The app header features the "FurLink" logo and "Pet CGM" text. The main heading is "Sign Up". The form includes an email input field with an envelope icon, a "Validation code" field, and a "Send" button. Below is a "Country" dropdown menu. Password requirements are listed: "At least 6 characters", "Uppercase letters", "Lowercase letters", and "At least 1 digit". There are input fields for "Enter password" and "Confirm password", each with an eye icon for toggling visibility. A radio button is selected for "Agree Use of Terms & Privacy Policy". A blue "Confirm" button is at the bottom, with a link "Already have an account? Sign-in here" below it.

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